

**TechNews: TN#534**

**Distribution:** Internal External

**Product:** CS Imaging 8 Client / Server

**Validity:** NA

**Date:** 9 pt Jan 18<sup>th</sup>, 2022

**Writer:** Service Engineering

**Document Reference:** TN#534-20220118-EN

**Contact:** Carestream Dental - 4, rue F. Pelloutier - Croissy-Beaubourg - 77435 Marne-la-Vallée Cedex 2 (France)

## CSI8\_SQL\_Tables\_Remover: Err\_S\_CSIS\_10206 (DPMS only) / Migration from scratch

**Tool: CSI8\_SQL\_Tables\_Remover.exe**

### Description:

When a DPMS launches CSI, the information sent are stored in the SQL database (Name, first name, birthdate, and patient folder). **If the patient information already exist**, CSI detects a problem and this error appears

**(Err\_S\_CSIS\_10206). The solution** is to remove all SQL records with the CSI8\_SQL\_Tables\_Remover.

To prevent this error from occurring again, it is necessary to **check the path used by the DPMS on all workstations**. This path must be strictly identical (case sensitive, no mapped drive).

“**CSI8\_SQL\_Tables\_Remover.exe**” can be used also to delete all records in the SQL database and do again the migration from scratch manually.

This tool is only compatible with **CSI8 Client/Server version (v3.x.x)**. SSMS installation is not needed.



### How to use the SQL Table Remover.exe?

Download the tool from InfoTech in [CSI8>Download>Patch/Hotfix](#)

On the server, launch CSI8\_SQL\_Tables\_Remover.exe and confirm the operation with “Y”.

Admin password can be requested.

**WARNING: Be careful, the tool will remove all SQL tables, no more patients will be displayed in CSI8. This tool is designed specially to fix the problems with a DPMS and to reset all SQL tables to do again the migration manually from scratch.**

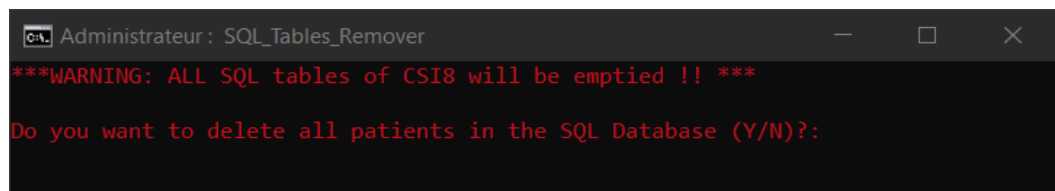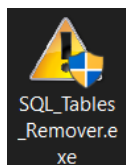

*Note: If all patients were deleted by mistake, open Monitor > CS Imaging Server > Service Tool and restore the last backup available. A backup was done when running the script.*